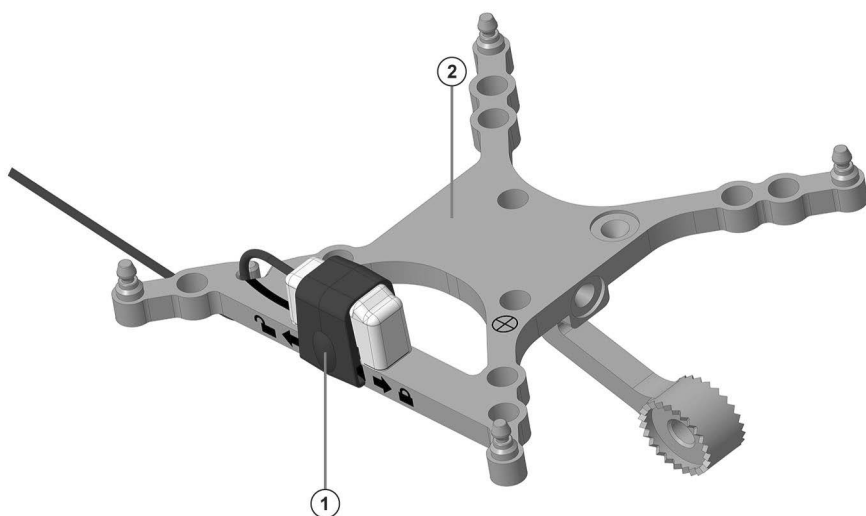


**EN**

EM Spine Reference Frame Tracker

(Component not shown to scale)

- 1 EM Spine Reference Frame Tracker (9736510)
- 2 Hybrid Reference Frame (9736511)

Device description

The electromagnetic (EM) Spine Reference Frame Tracker is a navigated device that consists of enclosed coils, a mounting fixture, and a communication cable. The EM Spine Reference Frame Tracker is a sterile, single-use device that is mounted to the Hybrid Reference Frame. The tracker coils are localized within the navigation system's emitter volume, generated by a low-energy magnetic field.

The EM Spine Reference Frame Tracker is only to be used with segmental tracking on the Stealth AXIS™ Surgical System. For instructions on using the Segmental Tracking Feature with the Stealth AXIS™ Surgical System, see the Spine clinical application instructions for use.

Intended use

The EM Spine Reference Frame Tracker is intended to enable location of patient anatomy during spinal procedures with Medtronic computer-assisted surgery systems.

Indications for use

The Stealth AXIS™ Surgical System is intended for precise positioning of surgical instruments and as an aid for precisely locating anatomical structures in open, minimally invasive, and percutaneous procedures. The Stealth AXIS™ Surgical System is indicated for medical conditions in which the use of stereotactic surgery may be appropriate, and where reference to a rigid anatomical structure, such as the skull or vertebra, can be identified relative to a CT or MR based model, fluoroscopy images, or digitized landmarks of the anatomy.

Contraindications

None.

Material composition

For materials of concern information such as REACH, CA Prop 65, or other product stewardship programs, go to www.medtronic.com/productstewardship.

Medtronic

Compatibility

The EM Spine Reference Frame Tracker is compatible with Stealth AXIS™ Surgical System running software version 1.0 and later.

System level manuals

This device can be used with Medtronic computer-assisted surgery systems. See the system level manuals and clinical application instructions for use for descriptions of patient groups, intended users, potential clinical benefits, side effects, and potential complications. This device is an accessory to your system.

See the system level manuals and the spine clinical application instructions for use for additional instructions, warnings, and precautions.

Content of this document

This document contains instructions for use on the EM Spine Reference Frame Tracker used during spinal procedures.

Warnings

Warning: Inspect sterile device packaging and expiration date before use. Do not use any pre-sterilized device whose sterile packaging appears to be compromised or whose expiration date has passed. Use of a compromised device could result in patient injury or death.

Warning: Prior to use, examine the device for damage and deterioration. Do not use any compromised device. Abandon use of any device damaged during the procedure.

Warning: Do not open sterile-barrier packages or containers until surgical use. At time of surgical use, inspect barrier for breach. If the sterile barrier is breached before surgical use, dispose of all components contained in the package.

Warning: Use personal protective equipment, as necessary, when handling, processing, and disposing of the device. Universal precautions for handling instruments after use should be observed by all hospital personnel according to OSHA Standard 29 CFR 1910.1030, Occupational Exposure to Blood borne Pathogens, or local regulations.

Warning: Do not cut the EM Spine Reference Frame Tracker cable before unplugging it from the system. Cutting the cable could pose an electrical safety impact to the patient.

Warning: Do not allow fluid to enter the device.

Warning: Post-registration movement of the EM Spine Reference Frame Tracker with respect to the patient reference frame or anatomy will result in inaccurate navigation, potentially resulting in patient harm. If you move the EM Spine Reference Frame Tracker at any time after registration, you must re-register.

Warning: Do not bump or reposition the EM Spine Reference Frame Tracker after registration. Movement of the EM Spine Reference Frame Tracker after registration might result in inaccurate navigation and could result in serious injury or death. If the device moves in relation to the patient anatomy at any time after assignment, you must detach and re-attach the EM Spine Reference Frame Tracker, and then verify accuracy by navigating on known anatomical landmarks. You must re-register if there is inaccurate navigation.

Warning: Do not replace the EM Spine Reference Frame Tracker while the Hybrid Reference Frame is attached to the patient. Replacing the EM Spine Reference Frame Tracker while the Hybrid Reference Frame is attached to the patient could move the anatomy and result in inaccurate navigation. If you must replace the EM Spine Reference Frame Tracker, remove the Hybrid Reference Frame from the patient, replace the EM Spine Reference Frame Tracker, and re-register the patient.

Warning: Do not move the EM Spine Reference Frame Tracker after patient registration. Movement of the EM Spine Reference Frame Tracker relative to its initial position on the Hybrid Reference Frame might result in inaccurate navigation and could result in serious injury or death. If you remove and replace the EM Spine Reference Frame Tracker at any time after registration, verify accuracy by navigating on known anatomical landmarks. You must re-register if there is inaccurate navigation. See the *Hybrid Reference Frame* instructions for use regarding reference frame placement considerations when using segmental tracking.

Warning: During navigation, visually confirm navigational accuracy frequently by localizing on known points. Check the displayed navigation information and verify that the system displays the name of the tracked tool. Abort use of the system if navigation seems inaccurate and steps to restore accuracy are unsuccessful.

Warning: This device is for single use only. After use, dispose of the device according to local and national regulations.

Precautions

Caution: Federal Law (U.S.A) restricts this device to sale by or on the order of a physician.

Caution: Avoid damage to the EM Spine Reference Frame Tracker including the cable. Do not use tracked devices if there are any signs of damage, deterioration, or deformation. If you damage the device during the procedure, discontinue use of the device immediately.


Use

Attach the EM Spine Reference Frame Tracker to the Hybrid Reference Frame before performing patient registration steps in the Spine clinical application. See the Spine clinical application instructions for use for information about patient registration.

Note: The EM Spine Reference Frame Tracker must be placed within the navigated field produced by the emitter. For additional information about restrictions on environments in which you can use the EM Spine Reference Frame Tracker, see the *Emitter Quick Reference Guide*.

1. Mount the EM Spine Reference Frame Tracker to the Hybrid Reference Frame.
 - a. Slide the top clip of the EM Spine Reference Frame Tracker back toward the cable connection before mounting onto the frame.
 - b. Locate the large and small posts on the crossbar of the Hybrid Reference Frame. Place the EM Spine Reference Frame Tracker over the two posts, ensuring the large post and the small post are aligned with the correct openings on the EM Spine Reference Frame Tracker.
 - c. Slide the clip toward the middle of the frame until it latches. The clip is properly latched when it snaps into place against the stop.
 - d. Thread the EM Spine Reference Frame Tracker cable through the notch in the frame for cable management.
 - e. Test the attachment of the EM Spine Reference Frame Tracker to the Hybrid Reference Frame by gently pushing and pulling on the sides. If the EM Spine Reference Frame Tracker does not move, the clip is properly latched.

Warning: Do not apply excessive force when attaching the EM Spine Reference Frame Tracker to the cable hook of the Hybrid Reference Frame. Applying excessive force could damage the cable.

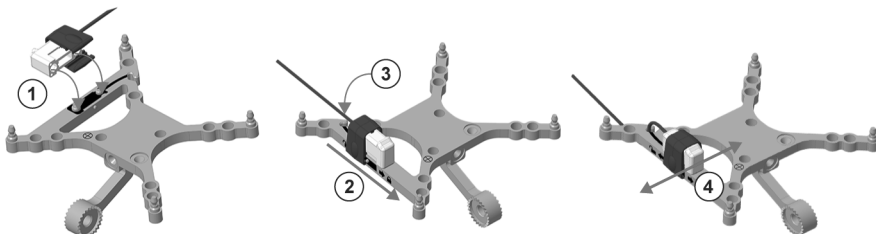

- 1 Mount the EM Spine Reference Frame Tracker to the posts
 - 2 Slide the clip into the locked position
 - 3 Thread cable through notch
 - 4 Apply gentle force to the sides of the EM Spine Reference Frame Tracker to check for secure attachment
 2. Connect the EM Spine Reference Frame Tracker cable to an open port on the EM instrument interface.

Note: Refer to the clinical application for information about how to assign the EM Spine Reference Frame Tracker to track desired vertebral levels.
 3. After the procedure, slide the clip toward the cable connection to release the latch and pull the EM Spine Reference Frame Tracker off the two posts to remove it from the frame.
 4. Dispose of the single use EM Spine Reference Frame Tracker.
- Note:** If you have questions about the device, contact Medtronic technical services.

Sterilization

Discard the device after use in accordance with the applicable regulations.

Warning: This device is for single use only. Do not reprocess, re-sterilize, or re-use this device. Attempts to reprocess or re-sterilize the device may be ineffective and may compromise its structural integrity. Any re-use creates a risk of contamination which could result in patient injury, illness, or death.

Note: This device is EO sterilized.

Disposal

Dispose of contaminated devices according to your facility's biohazard disposal procedure and local and national regulations.

Storage

Ensure the sterile packaging is dry before storing. Store in dry, clean conditions at ambient room temperature such that the package is not compromised.

Contact information

If a serious incident occurs in relation to the use of this device, report it to Medtronic Navigation. If a serious incident occurs in the European Union, also report it to the competent authority in the Member State where the incident occurred.

Assistance

If you have questions about these devices, contact your local Medtronic Navigation representative or call Medtronic Navigation, Inc. at 1 800 595 9709 or +1 720 890 3160.

Medtronic and Medtronic logo are trademarks of Medtronic. TM* Third-party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company.

This product may be protected by US Patents listed at: medtronic.com/patents.

Symbol definitions

The following symbols may appear on device packaging:



Caution: Federal law (U.S.A.) restricts this device to sale by or on the order of a physician.



Consult electronic instructions for use.



Type BF applied part.

R_x Only

©2026 Medtronic Navigation, Inc.

All Rights Reserved

M059985C001DOC1

Revision A

2026-02



Medtronic Navigation, Inc.

200 Medtronic Drive

Lafayette, CO 80026

USA

Main 1 720 890 3200

Fax 1 720 890 3500

www.medtronic.com

Technical support:

USA 1 800 595 9709

International +1 720 890 3160

rs.navtechsupport@medtronic.com

www.medtronic.com/manuals



M 05 99 85 C O O 1



R E V A