**EN**

EM Segment Tracker Driver

(Component not shown to scale)

1 EM Segment Tracker Driver (9736513)

Device description

The electromagnetic (EM) Segment Tracker Driver is an accessory to the Stealth AXIS™ Surgical System. The EM Segment Tracker Driver is a surgical instrument made of EM-compatible materials that is used to place segment trackers on patient anatomy during spinal procedures.

For instructions on using the Segmental Tracking Feature with the Stealth AXIS™ Surgical System, see the Spine clinical application instructions for use.

Intended use

The EM Segment Tracker Driver is intended to attach and detach the EM Segment Tracker during spine procedures.

Indications for use

The Stealth AXIS™ Surgical System is intended for precise positioning of surgical instruments and as an aid for precisely locating anatomical structures in open, minimally invasive, and percutaneous procedures. The Stealth AXIS™ Surgical System is indicated for medical conditions in which the use of stereotactic surgery may be appropriate, and where reference to a rigid anatomical structure, such as the skull or vertebra, can be identified relative to a CT or MR based model, fluoroscopy images, or digitized landmarks of the anatomy.

Contraindications

None.

Material composition

Material contained in product that can cause an allergic reaction: Nickel

Cobalt (CAS No. 7440-48-4) is classified as carcinogenic and toxic for reproduction substance of Category 1B. Affected products are labeled with Hazardous Material Symbol followed by 7440-48-4.

Stealth EM Segment Tracker Driver (9736513): Cobalt may be present at >0.1% weight/weight on direct patient contacting surfaces.

For additional materials of concern information such as REACH, CA Prop 65, or other product stewardship programs, go to www.medtronic.com/productstewardship.

Compatibility

The EM Segment Tracker Driver is compatible with EM Segment Trackers, which are compatible with the Stealth AXIS™ Surgical System.

System level manuals

This device can be used with Medtronic computer-assisted surgery systems. See the system level manuals and clinical application instructions for use for descriptions of patient groups, intended users, potential clinical benefits, side effects, and potential complications. This device is an accessory to your system.

See the system level manuals and the spine clinical application instructions for use for additional instructions, warnings, and precautions.

Medtronic

Content of this document

This document contains use instructions for the EM Segment Tracker Driver in segmental tracking procedures.

Warnings

- Warning:** Before first use and for reuse, clean and sterilize the device according to the instructions in this document. The device is designed for reuse and is supplied non-sterile. Use of a non-sterilized device could result in patient injury or death.
- Warning:** Use caution and exercise care to prevent injury during use. See the corresponding system level manuals and clinical application instructions for use for further handling instructions.
- Warning:** Do not open sterile-barrier packages or containers until surgical use. At the time of surgical use, inspect barrier for breach.
- Warning:** Abandon use of the device if the device sterility is compromised during the procedure. Use of a non-sterilized device could result in patient injury or death.
- Warning:** Do not hammer the segment tracker into the patient's anatomy. Hammering the tracker into bone could result in unstable fixation and can damage the tracker, and will affect navigational accuracy.


Precautions

Caution: Federal law (U.S.A.) restricts this device to sale by or on the order of a physician.

Use

Follow these steps to use the driver.

1. Connect the segment tracker into the EM instrument interface using the connector portion of the device. Ensure the EM instrument interface is connected to the workstation.
Note: If you want to navigate the tracker into place, complete patient registration before attaching the segment tracker to the patient.
2. Insert the segment tracker into the driver opening until the assembly is secure and fully seated.
3. Press the cable into the driver slot. Gently pull the cable taut to ensure the tracker is fully seated within the driver. Use the cable capture feature at the top of the driver to secure the cable in the driver.
Warning: Align the axis of the driver with the axis of the segment tracker screw when attaching the segment tracker to the patient anatomy. Misaligning the axes could damage the tracker and cause inaccurate navigation. Ensure that the tracker is fully seated into the driver and that there are no particulates inside of the driver. If the tracker is not properly seated in the driver, it could damage the tracker and cause inaccurate navigation.

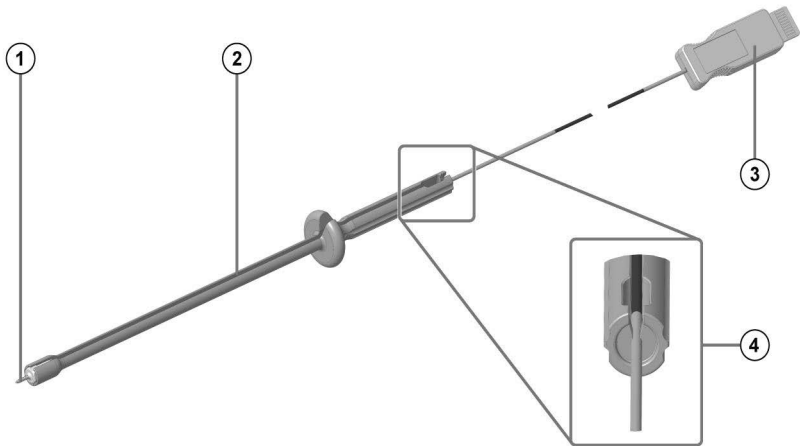


- | | |
|---|--|
| 1 | Segment tracker |
| 2 | Driver slot |
| 3 | Connector portion of the segment tracker |
| 4 | Cable capture feature |

4. Attach the segment tracker into the bone using the driver.

- a. To prevent cable wrapping, ensure that the hand turning the driver is underneath the cable and that there is enough cable in the sterile field to allow for sterile cable management throughout the procedure. Use the other hand to stabilize and guide the placement of the driver.
- b. To minimize the chances that the tracker will slip off the bone or that you will damage the tracker, align the driver perpendicular to the bone.
- c. Once the assembled tracker and driver are atop the desired tracker location, turn the driver clockwise while applying gentle downward force until the tracker is fully seated on bone.
- d. If you want to track multiple vertebrae with a single tracker, consider how you want to group levels when placing a tracker on a specific level. For instructions about how to group multiple vertebral levels to one segment tracker, see the clinical application.
- e. When attaching the EM Segment Tracker with navigation, see the *Emitter Quick Reference Guide* for emitter placement.


Warning: Do not apply excessive force when attaching the segment tracker to the patient. Applying excessive force could strip the bone and cause inaccurate navigation.

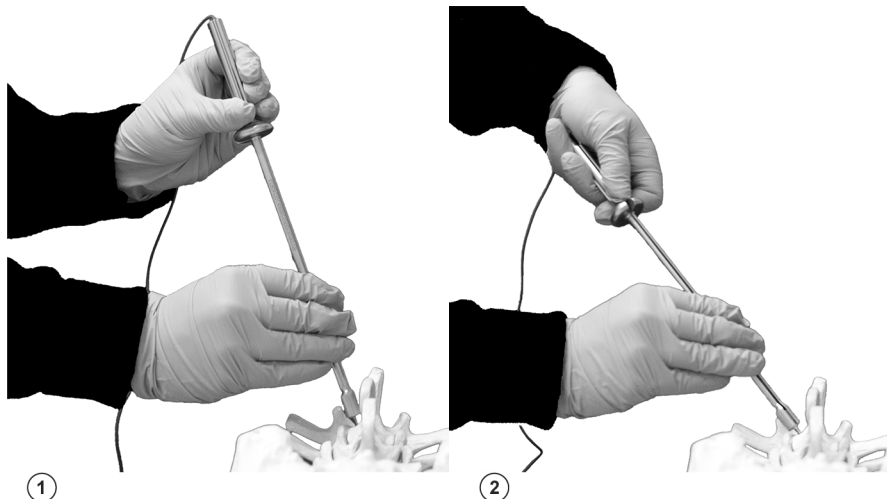

- 1 Hand position to gently secure tracker in bone with minimal pressure.
- 2 Hand position to apply gentle downward pressure if needed to gain purchase and progress screw into bone.

5. When the segment tracker is secured in bone, gently remove the cable from the driver slot. Gently pull up to remove the driver from the segment tracker.

Note: If you have questions about inserting the segment tracker into patient anatomy using the driver, contact Medtronic technical services.

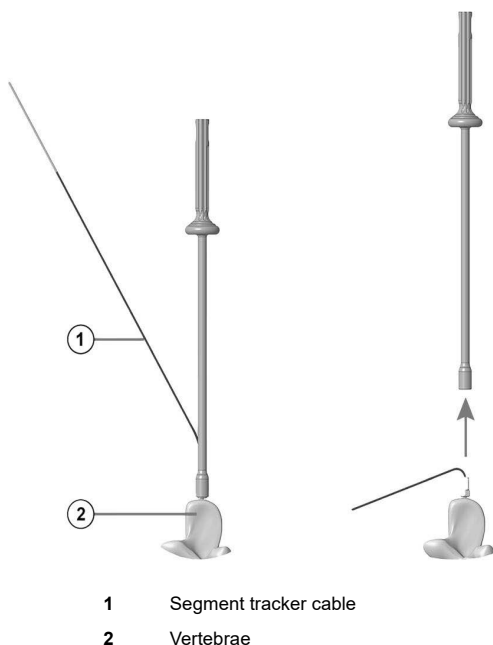

6. After the procedure, reattach the driver to the segment tracker and turn counterclockwise to remove the tracker from the bone.

Warning: Ensure that the tracker is fully seated into the driver and that there are no particulates inside of the driver. If the tracker is not properly seated in the driver, it could damage the tracker and cause inaccurate navigation or patient harm.

Note: Ensure the segment tracker is connected to the EM instrument interface until the segment tracker is removed from patient anatomy.

Troubleshooting

- If the segment tracker cable wraps around the driver during insertion, verify that the cable is secured in the cable capture feature on the end of the driver and that the cable is exiting the end of the driver with a clear path. You can also have another person hold the cable exiting the end of the driver while inserting the segment tracker.
- If you use excessive force while attaching the tracker in bone, you might strip the bone. Decrease downward force on the driver and use two or three fingers to rotate the driver. Do not use your entire hand to rotate the driver.
- If the segment tracker breaks, immediately remove the broken tracker and replace it with a new tracker.
- If you have trouble successfully attaching the tracker, try inserting the tracker into a different part of the bone.
- If it feels like the segment tracker is spinning freely within the bone, attach the tracker to a different part of the bone. Do not reuse holes to insert the segment tracker as this can lead to unstable fixation and inaccurate tracking.

Note: Different bone densities may require different amounts of effort to attach the tracker to bone.

Device processing - cleaning, inspection and sterilization

These instructions are intended to assist hospital personnel in safely handling and effectively processing devices supplied non-sterile by Medtronic. The specified devices must be cleaned, inspected, packaged, and sterilized prior to initial use and each subsequent use.

General warnings and precautions	Warning: Use CE-marked processing equipment that has been properly validated, calibrated, and maintained. Warning: To prevent cross-contamination, always handle, transport, clean, sterilize, and package devices contacting the central nervous system separately from other devices. Warning: Ensure that all processing steps are properly and completely executed before releasing the device for use. Warning: Do not use the device if there are any signs of damage, deterioration, or deformation. Discontinue use of the device immediately if it becomes damaged during the procedure.
Limitations on reprocessing	The useful life of this device depends on a variety of factors including method and duration of each use, as well as handling between uses. End of useful life is normally determined by wear and damage due to use. Careful inspection and functional testing prior to use is the best method for determining if the device is at the end of its useful life. See the "Maintenance, inspection, and testing" section below for additional instruction. Warning: Prior to use, examine the device components for damage, deterioration, deformation, and abuse. Do not attempt to use any device that appears to be bent or otherwise damaged. Abandon the use of any device damaged during the procedure. Refer to the "Maintenance, inspection, and testing" instructions in this document to determine if the device is suitable for use. If defects are found, do not use the product, and contact Medtronic immediately. Warning: Do not reuse or reprocess a device where prion contamination is suspected or confirmed. Incinerate the potentially contaminated device according to national regulations. Warning: At the end of its useful life, dispose of the device in accordance with national regulations.

Instructions

Initial treatment at the point-of-use	1. Processing begins at the point of use. Do not allow blood, debris, or bodily fluids to dry on the device. 2. Remove gross contamination from the device using non-linting or non-shedding wipes or cloths. 3. If the device cannot be reprocessed immediately, keep the device moist during transport by placing the device in a tray of utility water, if available, or by covering the device with a damp towel. 4. To prolong the life of the device, process immediately after use. Contain and transport the device in a closed, puncture proof container per OSHA or local regulations. Soiled devices should be transported separately from non-contaminated devices to avoid contamination of personnel and surroundings. Warning: Use personal protective equipment, as necessary, when handling, processing, and disposing the device. Universal precautions for handling devices after use should be observed by all hospital personnel according to OSHA Standard 29 CFR 1910.1030, Occupational Exposure to Bloodborne Pathogens, or local regulations. Warning: Devices should be cleaned within 30 minutes of use to limit fixation of contaminants.
Preparation before cleaning	1. Thoroughly rinse the device using cold running utility water (10-22°C [50-72°F]) for at least 30 seconds and until all visible soil is removed. While rinsing, actuate any movable components through their full range of motion (if applicable) and use a soft-bristled nylon brush (such as Key Surgical, part number N-2100) and clean non-linting wipe to aid in soil removal. Special attention should be given to any lumens, blind holes, mated surfaces, crevices, cavities and other features that may present a challenge to cleaning. 2. Carefully inspect the device to ensure that all visible soil has been removed. If any soil remains on the device, continue rinsing and brushing the device until it is visibly clean. Warning: Do not use metal brushes, scouring pads, or abrasive cleaners to clean the device as this may damage the device. Caution: Do not use processing agents that may cause corrosion (such as saline or cleaning agents that contain aldehyde, active chlorine, bromine, or iodide).

Automated Cleaning	1. Transfer the device to a washer-disinfector complying with ISO 15883. The device should be positioned such that all surfaces are easily accessible, device lumens, blind holes, and cavities are facing down, and hinges are open (to allow for draining). 2. Select the appropriate cycle parameters based on detergent type (i.e., enzymatic or alkaline cleaning cycle). 3. Upon completion of the automated wash cycle, carefully inspect all parts of the device for visible soil. If any soil is observed, repeat the cleaning process. Warning: Use the rigid container and instrument tray for sterilization only. Instruments must be cleaned separate from the container or tray. Warning: Do not place heavy instruments on top of delicate devices to prevent damage to the devices. Note: Cleaner and detergent manufacturers provide concentrated solutions, which must be diluted per the manufacturer's instructions for cleaning of medical devices.		
	Enzymatic cleaning cycle		
	Phase	Recirculation/ Soak Time	Water Temperature
	Prewash	2 minutes (02:00)	Cold utility water 10-22°C (50-72°F)
	Enzyme wash	2 minutes (02:00)	Hot utility water 43-55°C (109-131°F)
	Wash	2 minutes (02:00)	Hot utility water 66°C (151°F) (setpoint)
	Rinse	15 seconds (00:15)	Hot utility water 43-60°C (109-140°F)
	Thermal disinfection	1 minute (01:00)**	Hot critical water ≥ 90°C (194°F)
	Dry	30 minutes (30:00)	100°C (212°F) (setpoint)
	Note: Enzymatic cycle parameters were validated using STERIS Valsure™ Enzymatic Cleaner and Neutral Detergent. **Thermal disinfection has been validated for 1 minute at 90°C (A₀ ≥ 600) and for 5 minutes at 90°C (A₀ ≥ 3000).		
	Alkaline cleaning cycle		
	Phase	Recirculation/ Soak Time	Water Temperature
	Prewash	2 minutes (02:00)	Cold utility water 10-22°C (50-72°F)
	Alkaline wash	5 minutes (05:00)	Hot utility water 43°C (109°F) (setpoint)
	Rinse	1 minute (01:00)	Hot utility water 43-60°C (109-140°F)
	Thermal disinfection	1 minute (01:00)**	Hot critical water ≥ 90°C (194°F)
	Dry	30 minutes (30:00)	100°C (212°F)
	Note: Alkaline cycle parameters were validated using Dr. Weigert Neodisher™ Mediclean Forte. **Thermal disinfection has been validated for 1 minute at 90°C (A₀ ≥ 600) and for 5 minutes at 90°C (A₀ ≥ 3000).		

Drying	If a dry phase is not included in the washer-disinfector cycle: • Dry the device thoroughly inside and out to prevent rusting and malfunction. • Dry the device with a clean, soft, lint-free towel. • Special attention should be given to threads, hinges, crevices, cavities, and other areas where fluid can accumulate. • Open and close the device (if applicable) so that all areas are accessible. • Dry lumens, blind holes, and mated surfaces using instrument grade compressed air.																												
Maintenance, inspection, and testing	1. Visually inspect the device under good light conditions for completeness, damage, and excessive wear. Consideration should be given to the following: - All device components should be accounted for. - Laser etchings, engravings, and other markings should be legible. - No chips, cracks, or sharp-edged gouges should be present. - No bending, wearing, warping or pitting should be present. - No discoloration, rust, or corrosion should be present. 2. If any of the criteria listed above is not met, then the device is at the end of its useful life and should be disposed of according to national regulations. Warning: Do not process or use any device that has obvious damage or corrosion.																												
Packaging	The device should be packaged for sterilization by placement in a Medtronic-branded sterilization tray within a Medtronic-branded reusable rigid sterilization container (9736633, 9736632, 9736655, or 9736654) with disposable filters. Warning: Do not place additional devices inside the sterilization tray as this may impact sterilization efficacy and dry time. Warning: Do not wrap reusable rigid sterilization containers. The device may not be sterilized to an adequate sterility assurance level. Note: The user facility is responsible for using only accessories (such as chemical indicators, and biological indicators) that have been approved by the Food and Drug Administration for medical facilities inside the USA and its territories or conform to EN ISO 11607 for medical facilities outside the USA and its territories.																												
Sterilization	1. Load the packaged device into the steam sterilizer according to the sterilizer's operating manual. 2. Program the steam sterilizer to deliver one of the following validated pre-vacuum (dynamic air removal) steam (moist heat) sterilization cycles: <table><tr><th>Cycle Type</th><th>Temperature</th><th>Exposure Time</th><th>Minimum Dry time*</th></tr><tr><td colspan="4">For Healthcare Facilities in any geography</td></tr><tr><td>Pre-vacuum</td><td>132°C (270°F)</td><td>4 minutes (04:00)</td><td>30 minutes (30:00)</td></tr><tr><td>Pre-vacuum</td><td>135°C (275°F)</td><td>3 minutes (03:00)</td><td>30 minutes (30:00)</td></tr><tr><td colspan="4">For Healthcare Facilities outside the USA and its territories</td></tr><tr><td>Pre-vacuum</td><td>134°C (273°F)</td><td>3 minutes (03:00)</td><td>30 minutes (30:00)</td></tr><tr><td>Pre-vacuum</td><td>134°C (273°F)</td><td>18 minutes (18:00)</td><td>30 minutes (30:00)</td></tr></table> *Minimum dry times were validated using sterilizers with vacuum drying capabilities. Additional time may be required to dry the device under ambient atmospheric pressure or variable environmental conditions. Warning: Do not sterilize and supply for surgical use any device that is not visibly clean and free of particulates. If particulates are present, repeat reprocessing, starting with the Preparation before cleaning step. Warning: The device may not be sterilized to an adequate sterility assurance level without prior cleaning. Warning: Do not stack containers or trays during sterilization as this may limit steam penetration and prevent effective sterilization. Note: Steam for sterilization should be generated from water that has been treated to remove total dissolved solids, filtered to remove contaminants and water droplets, and supplied via piping without deadlegs or other stagnant zones where contamination might collect. Steam saturation should be greater than 97%.	Cycle Type	Temperature	Exposure Time	Minimum Dry time*	For Healthcare Facilities in any geography				Pre-vacuum	132°C (270°F)	4 minutes (04:00)	30 minutes (30:00)	Pre-vacuum	135°C (275°F)	3 minutes (03:00)	30 minutes (30:00)	For Healthcare Facilities outside the USA and its territories				Pre-vacuum	134°C (273°F)	3 minutes (03:00)	30 minutes (30:00)	Pre-vacuum	134°C (273°F)	18 minutes (18:00)	30 minutes (30:00)
Cycle Type	Temperature	Exposure Time	Minimum Dry time*																										
For Healthcare Facilities in any geography																													
Pre-vacuum	132°C (270°F)	4 minutes (04:00)	30 minutes (30:00)																										
Pre-vacuum	135°C (275°F)	3 minutes (03:00)	30 minutes (30:00)																										
For Healthcare Facilities outside the USA and its territories																													
Pre-vacuum	134°C (273°F)	3 minutes (03:00)	30 minutes (30:00)																										
Pre-vacuum	134°C (273°F)	18 minutes (18:00)	30 minutes (30:00)																										
Storage	1. The device should remain in its sterilization packaging and be stored in a dry and clean location under ambient temperature conditions. 2. Consideration should be given to the sterile barrier system, storage manner, and environmental conditions when determining the maximum amount of time permitted between processing and use.																												

Additional information	• Definitions of and performance qualification levels for utility water and critical water are described in ANSI/AAMI ST108. • In the event of an alleged product failure, it is important to return the device to Medtronic for investigation; however, Medtronic is only able to accept devices that are visibly clean and free of contaminants.
------------------------	---

Medtronic has validated the processes provided in these instructions in accordance with ANSI/AAMI ST98, ISO 17665, AAMI TIR12 and ISO 17664 and demonstrated that they are capable of yielding product that meets quality, safety, and efficacy requirements. It remains the responsibility of the healthcare facility to ensure that processing, as performed using the equipment, materials, and operators at their facility, consistently achieves the desired result. This requires verification or validation or both, and routine monitoring of the process. Alternative methods of processing outside the scope of these instructions may be suitable; however, they must be validated by the end user. In the event of conflicting national cleaning or sterilization requirements, those shall prevail over Medtronic's recommendations.

Disposal

Dispose of contaminated devices according to your facility's biohazard disposal procedure and local and national regulations.

Contact information

If a serious incident occurs in relation to the use of this device, report it to Medtronic Navigation. If a serious incident occurs in the European Union, also report it to the competent authority in the Member State where the incident occurred.

Assistance

If you have questions about this device, contact your local Medtronic Navigation representative or call Medtronic Navigation, Inc. at 1 800 595 9709 or +1 720 890 3160.

Medtronic and Medtronic logo are trademarks of Medtronic. TM* Third-party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company.

This product may be protected by US Patents listed at: medtronic.com/patents.

Symbol definitions

The following symbols may appear on device packaging:

Rx Only

Caution: Federal law (U.S.A.) restricts this device to sale by or on the order of a physician.



Consult electronic instructions for use.

R_x Only

©2026 Medtronic Navigation, Inc.

All Rights Reserved

M060004C001DOC1

Revision A

2026-02



Medtronic Navigation, Inc.

200 Medtronic Drive

Lafayette, CO 80026

USA

Main 1 720 890 3200

Fax 1 720 890 3500

www.medtronic.com

Technical support:

USA 1 800 595 9709

International +1 720 890 3160

rs.navtechsupport@medtronic.com

www.medtronic.com/manuals



M 06 00 04 C 0 0 1



R E V A